


Procedure for Initial Adjustment

Button	Step	Action	Button	Step	Action
	1	If the pump is not already running, please turn on the pump.		6	Using the left-side + / - arrows, adjust to the desired percentage.
	2	If the "Lock" button is activated, press and hold it for 5 seconds to deactivate it.		7	Press "AP" once to save.
	3	Press the "Setup" button 2 times to display "CUSTOM" on the screen.		8	Using the left-side + / - arrows, adjust "AP" to 5 min.
	4	Using the right-side Firm / Soft arrows, increase or decrease the pressure according to the startup adjustment table.		9	To secure the settings, press "Lock".
	5	Press and hold "AP" for 5 seconds.	<i>The "LAL" (Low Air Loss) and "Max" buttons will remain accessible despite the lock symbol being active</i>		

If you do not know the user's weight, set to 28-16, 60% AP at 5 min (SEE PUMP SCREEN BELOW and procedure).

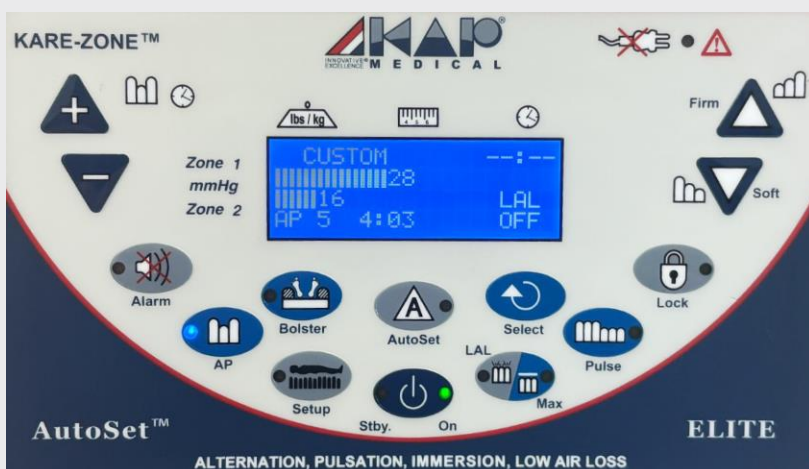

OTHER SETTINGS



Screen Lock

Press and hold the lock icon for 5 seconds. Only the "LAL / Max" button remains accessible in locked mode.



LAL (Low Air Loss)

Activate in case of excessive sweating or moisture



Max

Activate during care, transfers, and mobilization (automatically returns to alternating AP after 30 minutes for safety reasons)



Pulse

Wave effect, gentle stimulation through pulsation.

This mattress is equipped with the Seat Inflate (raised backrest) mode, an automatic inflation adjustment system designed to optimize support when the backrest is raised.



Startup Adjustments

K2M, K2MP, K2ME SURFACES | "CUSTOM" MODE

If you do not know the user's weight, set to :

	28	60%	28-16	5 min
Weight kg / lb	mmHg	%AP	mmHg zone 1 and mmHg zone 2	AP (duration)
23 kg / 50 lb	18	75	18 - 13	5 min
34 kg / 75 lb	20	75	20 - 15	5 min
45 kg / 100 lb	22	60	22 - 13	5 min
57 kg / 125 lb	23	50	23 - 11	5 min
68 kg / 150 lb	24	50	24 - 12	5 min
79 kg / 175 lb	25	50	25 - 12	5 min
91 kg / 200 lb	26	50	26 - 13	5 min
102 kg / 225 lb	27	50	27 - 13	5 min
113 kg / 250 lb	28	60	28 - 16	5 min
125 kg / 275 lb	29	60	29 - 17	5 min
136 kg / 300 lb	30	75	30 - 22	5 min
147 kg / 325 lb	32	75	32 - 24	5 min
159 kg / 350 lb	32	75	32 - 24	5 min

For any questions regarding the other settings, please do not hesitate to contact us.

