


LEIKA Fall guard system

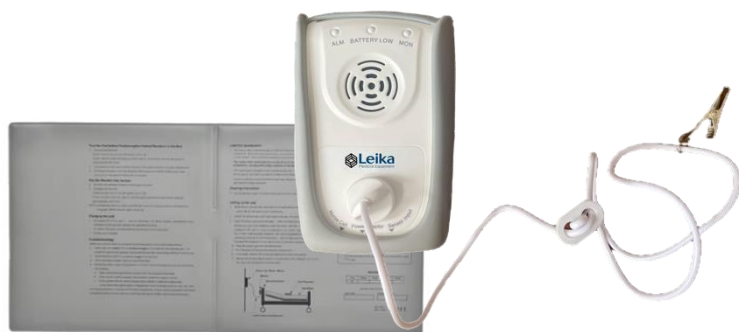

The LEIKA Fall guard system is a cord and voice patient safety monitor designed to detect resident movement and alert caregivers in real time. Compatible with both magnetic pull-cord and sensor cushion configurations, it combines reliable alarm technology with personalized voice message capability for a more reassuring patient experience.

The device features a customizable audio system that allows a familiar caregiver or family member's voice to reassure the patient and indicate that help is on the way. Automatic reset, nurse call compatibility, and a durable housing make it ideal for acute care, long-term care, and rehabilitation environments.

Main Features

- Compatible with magnetic pull-cord or sensor cushion
- Alarm triggered on movement: removal or magnet release
- Customizable voice message recording (up to 10 seconds)
- Nurse call system compatible
- Automatic reset
- Delay timer to reduce false alerts on brief movement
- HOLD function: repositioning window of 300 seconds
- 8 selectable alarm tones + 8 melodies from room/caregiver identification
- 3 alarm output modes: tone only, voice only, voice then tone
- Low battery notification with visual and audible alert
- Sensor loss notification
- No on/off switch – always ready
- Traction relief system protects sensor cushions from damage
- Thick, durable plastic housing
- Operates on 9V battery or optional AC power adapter
- 3 volume settings: Silent, 90 dB (Low), 108 dB (High) at 30 cm
- Mountable on wheelchair or bed
- Latex-free

Description

ALARM INDICATORS

- **ALM:** Flashing red light – alarm triggered
- **BATTERY LOW:** Flashing red +

Description

VOICE MESSAGE RECORDING

- Up to 10 seconds of custom voice recording
- Re-recordable at any time
- Familiar voice reassures patient while alerting staff

VOICE MESSAGE RECORDING

- Up to 10 seconds of custom voice recording
- Re-recordable at any time
- Familiar voice reassures patient while alerting staff

NURSE CALL INTEGRATION

- Compatible with Continuous Closure and One Shot nurse call systems
- Silent alarm mode available to reduce room noise

POWER

- 9V alkaline battery (6LR61)
- Optional AC power adapter

Other Information

MAINTENANCE

- Clean with damp cloth and mild detergent – do not immerse
- Avoid caustic cleaning products
- Inspect all components regularly; replace worn parts immediately
- Store in original box away from moisture and heat

